


Gastrostomy Tube Introduction Set

Non-endoscopic placement technique

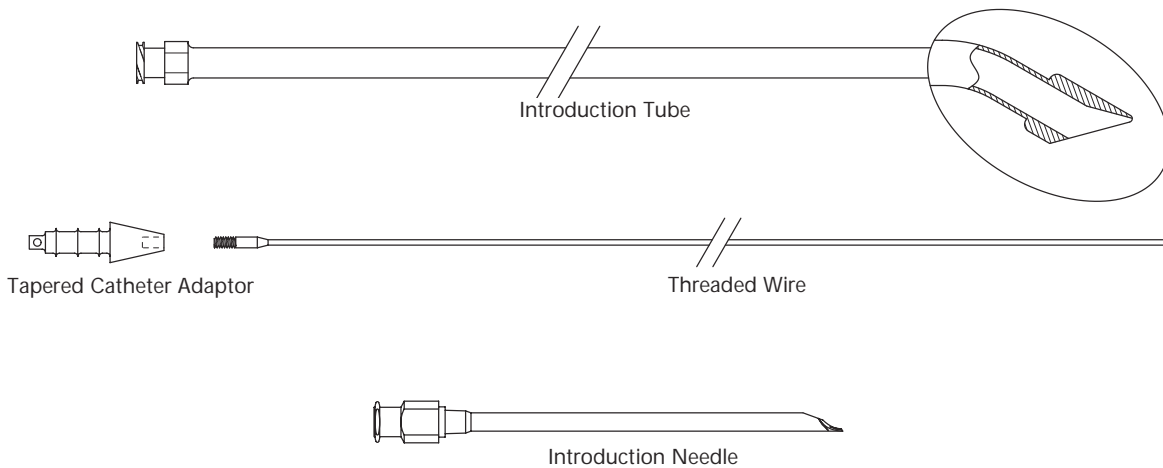

DESCRIPTION

- A reusable introduction set that enables the non-endoscopic placement of gastrostomy tubes into the stomach.
- Designed to accept either 16 or 20 French Cook® gastrostomy tubes.
- Intended for VETERINARY USE ONLY.

ORDER No.

Remarks

V-GTIS-1000

Set includes a 5.0mm O.D., 45cm long stainless steel stomach introduction tube with flared distal end, a 10 gauge, 5.0cm long percutaneous needle, a 70cm long 1.0mm stainless steel wire with threaded end and a tapered catheter adaptor. Recommended for animals up to 10kg (22lbs).

V-GTIS-2000

Set includes a 9.5mm O.D., 85cm long introduction tube with flared distal end, a 10 gauge, 5.0cm long percutaneous needle, a 110cm long 1.0mm stainless steel wire with threaded end and a tapered catheter adaptor. Recommended for animals over 10kg (22lbs).

Please note: Components can be ordered individually. Please contact your local Cook distributor for further details.

Gastrostomy/Oesophagostomy tubes...

Cook Veterinary Products' Gastrostomy/Oesophagostomy tubes are manufactured from silicone, and offer a number of advantages over similar tubes made from rubber or plastic:

- Silicone allows these tubes to be autoclaved and reused.
- Unlike other materials, silicone is resistant to degradation by stomach acid.
- Silicone is biologically inert, soft and flexible thus increasing patient comfort.