


# Flexi-Tip® Urethral Catheter

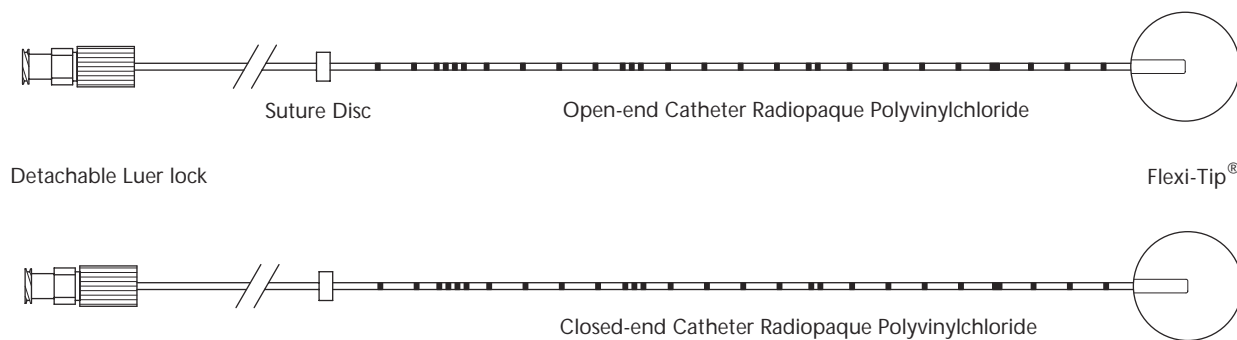

## DESCRIPTION

- A PVC catheter with a softer flexible tip for atraumatic cannulation of the urethra in dogs and other small animals.
- Available in closed or open ended versions.
- Distance centimetre markings on the catheter assist in accurate placement.
- Detachable Luer lock fitting allows the catheter to be trimmed to any length.
- Silicone movable disc on the catheter allows suturing to the prepuce to assist in catheter retention.
- Supplied sterile in peel-open packages.
- Intended for ONE-TIME USE.

ORDER No.                      Remarks

### CLOSED-END

V-021105                      5.0 French, 70cm closed-end Flexi-Tip® urethral catheter. Removable connector supplied.

### OPEN-END

V-021305                      5.0 French, 70cm open-end Flexi-Tip® urethral catheter. Removable connector supplied.

## The Solution to Difficult Introductions

Catheter introduction in dogs can be difficult due to the presence of strictures, uroliths, an enlarged prostate and the curve of the urethra as it passes around the ischial arch of the pelvis.

Flexi-Tip® urethral catheters have a flexible distal tip that allows them to negotiate the urinary tract in an atraumatic manner.

The catheter has 1cm markings to confirm how far the catheter has been introduced or where strictures or uroliths are positioned. A movable plastic Luer lock hub allow the catheter to be trimmed to a required length and the hub refitted.

### POSSIBLE APPLICATIONS

- Collection of urine samples for diagnostic purposes.
- Relieve urinary retention.
- Bladder lavage.