


Central Venous Catheters

Central Venous Catheter Set

with Peel-Away® Sheath Introduction Technique

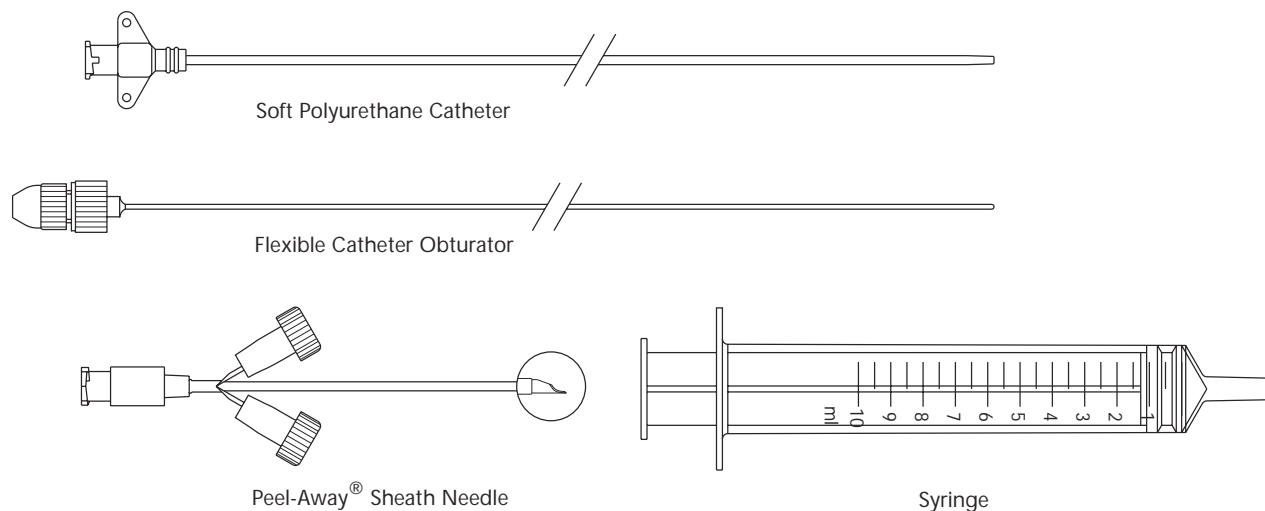

DESCRIPTION

- Used for venous pressure monitoring, blood sampling and administration of drugs and fluids.
- Central venous catheter introduced using a Peel-Away® sheathed needle.
- The Peel-Away® sheath preloaded on the access needle allows percutaneous introduction of the venous catheter, following direct venipuncture and removal of the needle cannula.
- After positioning of the catheter, the sheath is peeled away.
- Intended for ONE TIME USE.
- Intended for VETERINARY USE ONLY.


ORDER No.

Remarks

V-PPC-5.0U-15	Set consists of 5 French polyurethane catheter 15cm long, matched flexible obturator, 15 gauge, 7.0cm long Peel-Away® sheath needle and 10ml syringe.
V-PPC-4.0U-10	Set consists of a 4 French polyurethane catheter 10cm long, matched flexible obturator, 17 gauge, 7.0cm long Peel-Away® sheath needle and 10ml syringe.
V-PPC-3.0P-8	Set consists of 3 French polyethylene catheter 8cm long, a matched flexible obturator, 19 gauge, 7.0cm long Peel-Away® sheath needle and 10ml syringe.

CATHETER SET WITH PROXIMAL EXTENSION

Please note: Set includes flexible obturator, Peel-Away® sheath needle and syringe as show above.

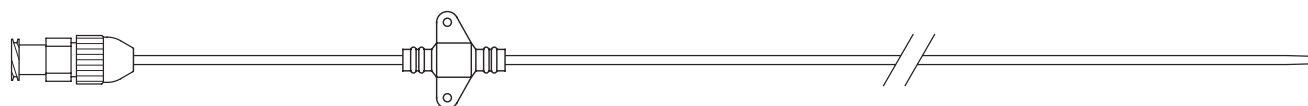

ORDER No.

Remarks

V-PPC-5.0U-20-PE5.0	Set consists of 5 French polyurethane catheter, 20cm long with 5.0cm proximal extension, matched flexible obturator, 15 gauge 7.0cm long Peel-Away® sheath needle and 10ml syringe.
---------------------	-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------