



Bone Biopsy Sets

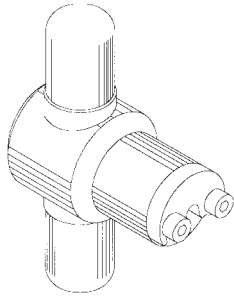

Biopsy Cutting Handle

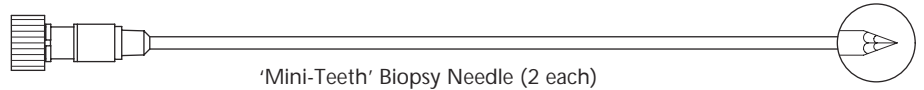

'Mini-Teeth' Biopsy Needle (2 each)

DESCRIPTION

- Used for biopsy of superficial and deep bone or bone marrow.
- The fine introducing needle guide enables safe positioning of guiding cannula through which multiple samples can be obtained with the 'mini-teeth' biopsy needles.
- The handle enables sufficient force to be placed on the biopsy needles for controlled biopsy of a variety of bone types.
- Supplied sterile in peel-open packages.
- Biopsy needles intended for ONE TIME USE.
- Intended for VETERINARY USE ONLY.




ORDER No.

Remarks

V-GBBS-200-6.0

Set consists of items shown above. Two 'mini-teeth' biopsy needles 14 gauge, 6.0cm long supplied.

Disposable Cytology Brushes

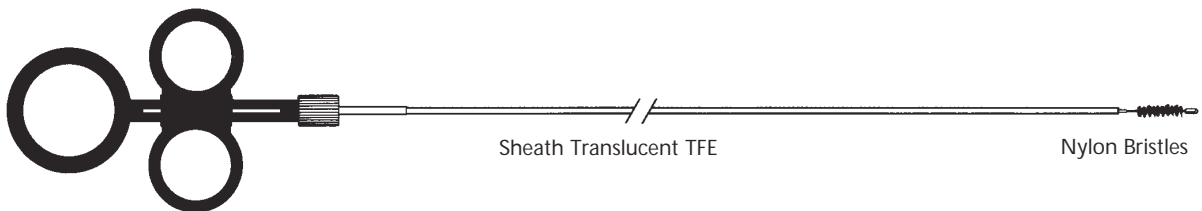

Three-Ring Handle Polycarbonate

Sheath Translucent TFE

Nylon Bristles

DESCRIPTION

- Used to collect cells from the bronchi or gastrointestinal tract for cytology.
- Ring handle enables the brush to be retracted into the sheath to prevent contamination of the sample.
- Supplied sterile in peel-open packages.
- Intended for ONE TIME USE.
- Intended for VETERINARY USE ONLY.



ORDER No.

Remarks

V-ECB-5-180-3-S

5 French, 180cm long cytology brush with 3.0mm brush diameter.